**Oral-B®**




**CrossAction®  
POWER MAX**

**B 1011**

**Type 3733  
Modèle 3733  
Modelo 3733**

99698216/XI-10

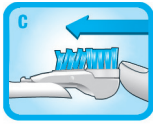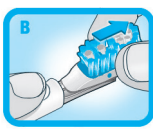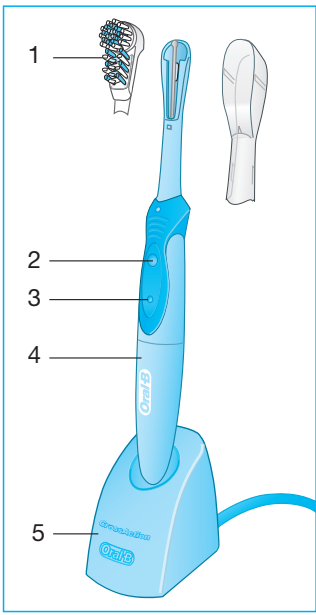Thank you for purchasing an Oral-B CrossAction PowerMAX. If you have any questions, please call:

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1-800-566-7252

# IMPORTANT SAFEGUARDS

When using electrical products, especially when children and persons with special needs are present, basic safety precautions should always be followed.

Read all instructions before using.

## DANGER

In order to reduce the risk of electrocution:

1. Do not place or store the charger where it can fall or be pulled into a tub or sink.
2. Do not place the charger in water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.
4. Do not use while bathing.

---

## WARNING

To reduce the risk of burns, electrocution, fire or injury:

1. Close supervision is necessary when this product is used by, on, or near children or invalids.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Keep cord away from heated surfaces.
4. Never drop the appliance or insert any object into any opening of the appliance.
5. Never operate if any part of this product has been dropped or damaged, if the cord or plug is damaged or the charger has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.

## IMPORTANT

1. Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Center. A damaged or non-functioning unit should no longer be used.
2. If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
3. Do not touch the power plug with wet hands. This can cause an electric shock.

4. When unplugging, always hold the power plug instead of the cord.
5. Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Service center. Do not disassemble the product except when disposing of the battery.
6. When taking out the battery for disposal, use caution not short the positive (+) and negative (-) terminals.
7. This product is not intended for use by children under age three.

## CAUTION

If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

## SAVE THESE INSTRUCTIONS

---

### Description

- 1 Brush head
- 2 On switch
- 3 Off switch
- 4 Handle
- 5 Charging unit

### Specifications

Power supply: see bottom of the charging unit.

Handle voltage: 1.2 V

### Connecting and charging

- Plug the charging unit into an electrical outlet. Place the handle on the charging unit. The handle must be switched off.
- A full charge takes 16 hours and will provide an operation time of approximately 60 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible.
- However for environmental reasons, Oral-B recommends to unplug the charging unit until the next charge is required.

### Battery maintenance

To maintain the maximum capacity of the rechargeable battery, you must discharge it at least once every 6 months. To do this, unplug the charging unit and allow the battery to run down with normal use.

### Using the toothbrush

- To avoid splashing, guide the brush head to the teeth, before switching the appliance on.
- Brush the gums as well as the teeth, first the outsides,

then the insides and finally the chewing surfaces.

- Brush for at least 2 minutes.

## Brush head replacement

1. **To remove:** Place QuickRelease™ tool\* on neck of toothbrush as shown and press (A), then slide brush head off toothbrush (B).
2. **To replace:** Lay new replacement brush head parallel to neck of toothbrush (C), then slide head toward handle until lock clicks into place (D). Check that brush head is firmly attached. Do not insert brush head at angle (E).

Change your brush head when the blue Indicator® bristles fade halfway or every 3 months, as recommended by dentists or dental hygienists.

- \* included in CrossAction® Power replacement brush.

## Environmental notice

### Nickel-Cadmium Battery

SAVE THIS INFORMATION  
PRODUCT CONTAINS A  
NICKEL-CADMIUM  
RECHARGEABLE BATTERY.  
MUST BE RECYCLED OR  
DISPOSED OF PROPERLY.

May not be disposed of in municipal waste.  
Recycle or dispose of in an environmentally sound manner.

Ni-Cd

Subject to change without notice.

**Limited warranty:** In case of defect in materials or workmanship, the appliance will be replaced within 2 years of date of purchase. If the appliance exhibits such a defect call 1-888-287-9050 for address to return appliance. We will send you a replacement appliance free of charge. Proof of date of original purchase is required. This warranty does not cover use and wear of batteries and brush head, nor does it cover damage to finishes, normal wear, or damage resulting from accident, misuse, dirt, tampering, or unreasonable use.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF FITNESS FOR ANY PARTICULAR PURPOSE, ARE EXPRESSLY DISCLAIMED. IN NO EVENT WILL MANUFACTURER BE RESPONSIBLE FOR CONSEQUENTIAL; INCIDENTAL OR SPECIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states/provinces do not allow the exclusion or limitation

of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state/ province to province.

**To obtain service in the USA**

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Braun Service Center closest to you.

To order original Braun replacement parts and accessories (including replacement brush heads) by phone:  
Please call 1-800-566-7252.

To speak to a Braun Consumer Service Representative:  
Please call 1-800-566-7252.

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