
PULSONIC

SLIM

powered by BRAUN

S 15.523.2

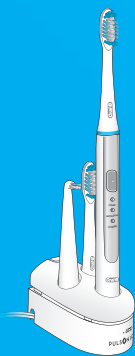English	3
Español	8
Français	12








Internet :

www.oralb.com
www.braun.com
www.service.braun.com
www.oralb.ca

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99434617/V-10
GB/E/F

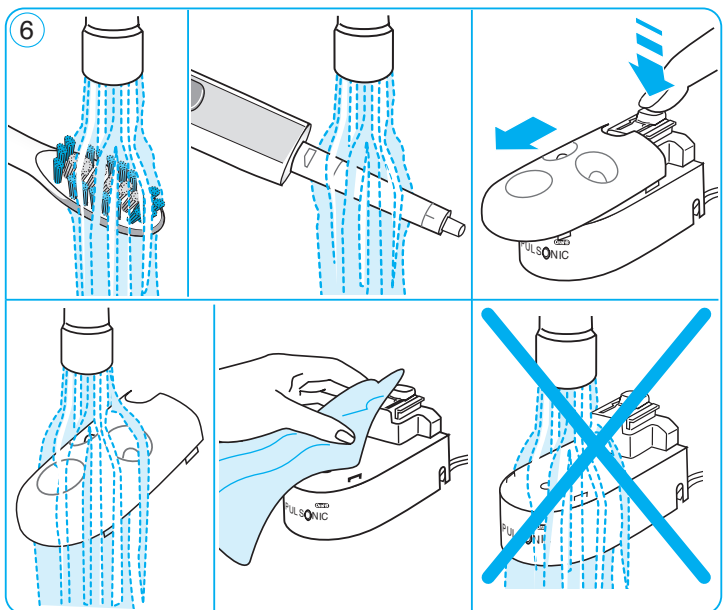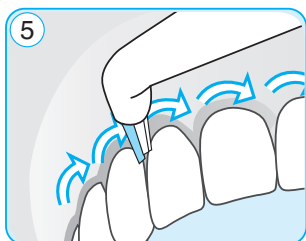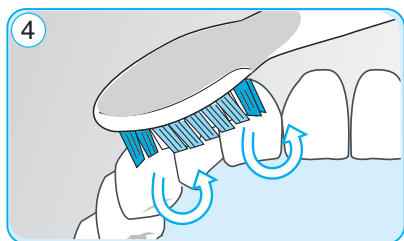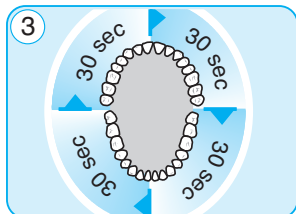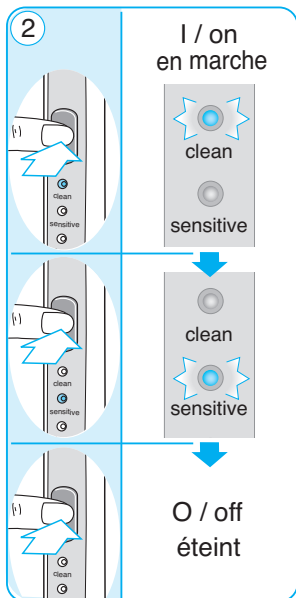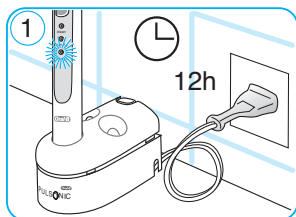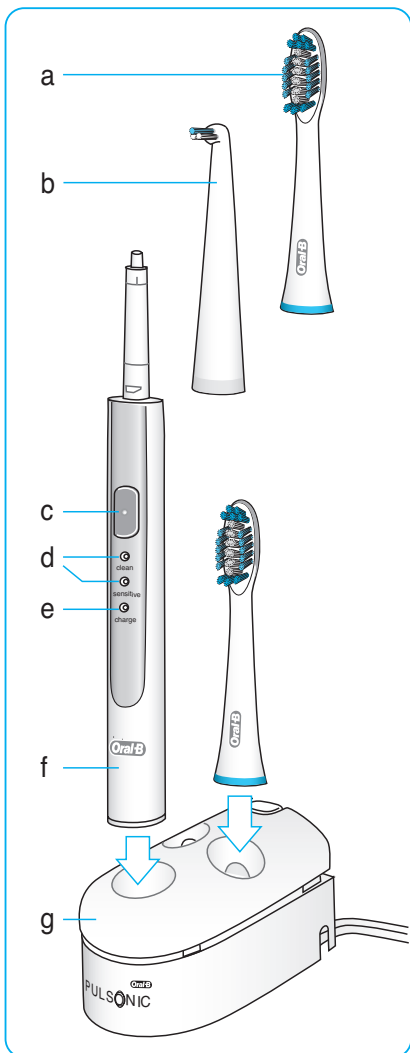IMPORTANT SAFEGUARDS

Oral-B Pulsonic Slim has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective. When using electrical products, however, basic safety precautions should always be followed.

Read all instructions before using.

DANGER

In order to reduce the risk of electrocution:

1. Do not place or store the charger where it can fall or be pulled into a tub or sink.
2. Do not place the charger in water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.
4. Please do not use while bathing.

WARNING

To reduce the risk of burns, electrocution, fire or injury:

1. Close supervision is necessary when this product is either used by or near children or invalids.
2. Use this product only for its intended use as described in this manual. Do not use attachments or chargers not recommended by the manufacturer.
3. Keep cord away from heated surfaces.
4. Never drop the appliance or insert any object into any opening of the appliance.
5. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged or if it has been dropped into water. If any of these occur, please contact a service centre and return the product for examination and repair.

IMPORTANT

1. Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral B Braun Service Center. A damaged or non-functioning unit should no longer be used.
2. Do not touch the power plug with wet hands. This can cause electric shock.
3. When unplugging, always hold the power plug instead of the cord.
4. Do not put metal objects such as coins and clips on the charger. This may generate heat damaging the unit and heating the metal item presenting a potential burn hazard.
5. Do not modify the product. Do not disassemble or repair. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Braun Service centre. Do not disassemble the product except when disposing of the battery.
6. When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
7. This product is not intended for use by children under age three.

CAUTION

If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

SAVE THESE INSTRUCTIONS

Description

- a Oral-B Pulsonic brush head
- b Oral-B Precision Tip brush head
- c On/off button
- d Mode indicators
- e Charge indicator
- f Handle
- g Charging unit with handle stand and release button

Connecting and charging

Your toothbrush is electrically safe and designed for use in the bathroom.

- Plug the charging unit (g) into an electrical outlet and place the handle (f) on the charging unit (see picture ①, page 2).
- A full charge takes at least 12 hours and enables approximately one week of regular brushing (twice a day, 2 minutes).
- For everyday use, the toothbrush handle can be stored on the plugged-in charging unit to maintain it at full power. Overcharging is impossible. However for environmental reasons, Oral-B recommends to unplug the charging unit until the next recharge is required.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

Brushing

Brushing modes

Oral-B Pulsonic Slim offers different brushing modes for varying oral care needs:

«Clean» – Gentle, effective mouth cleaning.

«Sensitive» – Gentle cleaning for sensitive areas.

When pressing the on/off button, your toothbrush automatically starts in the «Clean» mode. To shift from the «Clean» to the «Sensitive» mode, press the on/off button once again. To turn your toothbrush off, press the on/off button (see picture ②, page 2).

Brushing technique

- Slide the brush head on the handle. Wet brush head and apply toothpaste; you can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching on your toothbrush.
- Place the brush head on the outside surface of your upper teeth. The toothbrush bristles should be placed against the teeth at a slight angle towards the gumline.
- Turn your toothbrush on and start brushing in a slightly circular motion (see picture ④, page 2). Apply light pressure during brushing. After a few seconds, guide the bristles to the next section. Brush inside, outside and chewing surfaces of your teeth with the same motion throughout the mouth for the recommended 2 minutes.
- To avoid splashing, turn your toothbrush off before removing from the mouth by pressing the on/off button.

You may experience a slight tickling or tingling sensation when you first use your Oral-B Pulsonic Slim. This will subside as your teeth and gums will become accustomed to your Oral-B Pulsonic Slim cleaning experience. If you experience this initial tingling, you may want to brush in the «Sensitive» mode first.

In the first few days of using your toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or your dental hygienist.

Brush heads (depending on model)

Oral-B Pulsonic brush head

The Oral-B Pulsonic brush head (a) is recommended for everyday deep cleaning. Its bristles provide a thorough cleaning between your teeth, tooth surfaces and along the gumline. The blue indicator bristles help you to monitor brush head replacement need. With proper brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months, signalling it is time to replace your brush head.

Oral-B Precision Tip

The Oral-B Precision Tip (b) is designed to clean between teeth as well as around bridges, crowns and implants. To use, position it at a 45-degree angle to gumline and slowly guide it along the gumline from one interdental space to the next (see picture ⑤, page 2).

Timer

A short stuttering sound at 30-second intervals (one break) reminds you to brush all four quadrants of your mouth equally (see picture ③, page 2). A long stuttering sound (two breaks) indicates when you reach the professionally recommended 2-minute mark.

Cleaning


Switch off the handle and remove the brush head. Rinse both brush head and handle thoroughly under running water. Wipe them dry before you put the handle back on the charging unit. From time to time you can also clean the charging unit with a damp cloth. For thorough cleaning you may remove the handle stand by pushing the release button on top of the charging unit. The charging unit should never be placed in water (see picture ⑥, page 2).

Environmental notice

This product is provided with a nickelhydride battery, which does not contain any harmful heavy metals. Nevertheless, in the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at approved recycling or disposing locations in accordance with local regulations.

Battery removal at the end of the product's useful life

Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

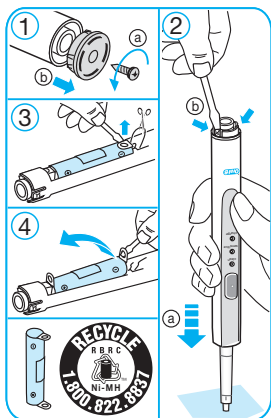Caution! Opening the handle will destroy the appliance and invalidate the warranty.

Subject to change without notice.

FOR USA ONLY

Limited 2-Year Warranty

Oral-B Braun warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Oral-B Braun will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Oral-B Braun Service Centre. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Centre is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Oral-B Braun.

This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL ORAL-B BRAUN BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

How to obtain service in the USA

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Braun Service Centre closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Braun Service Centre of your choice or,
- B. Ship the product to the authorized Oral-B Braun Service Centre of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Oral-B Braun replacement parts and accessories:

Contact an authorized Oral-B Braun Service Centre. Please call 1-800-566-7252 to be referred to the authorized Oral-B Braun Service Centre closest to you.

To order original Oral-B Braun replacement parts and accessories (including replacement brush heads) by phone:

Please call 1-800-566-7252.

To speak to an Oral-B Braun Consumer Service Representative:

Please call 1-800-566-7252.

FOR CANADA ONLY

Limited 2-Year Warranty

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service center, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages.

The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

How to obtain service in Canada

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Braun Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Braun Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Braun Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Oral-B Braun replacement parts & accessories:

Contact an authorized Oral-B Braun Service Center. Please call 1-800-566-7252 to be referred to the authorized Oral-B Braun Service Center closest to you.

To speak to an Oral-B Braun Service Representative:

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