




Oral-B®

VITALITY

powered by **BRAUN**

Type/Modèle/Modelo 3709

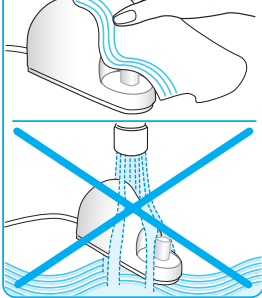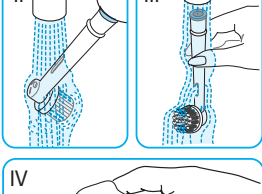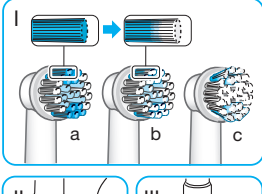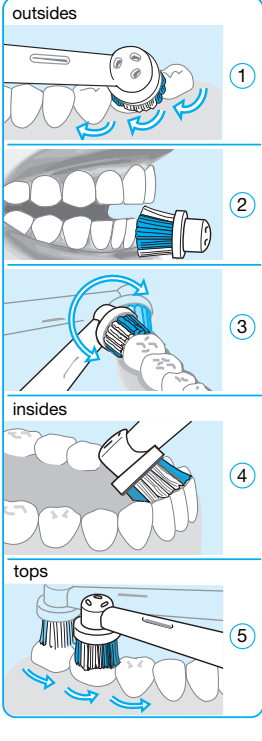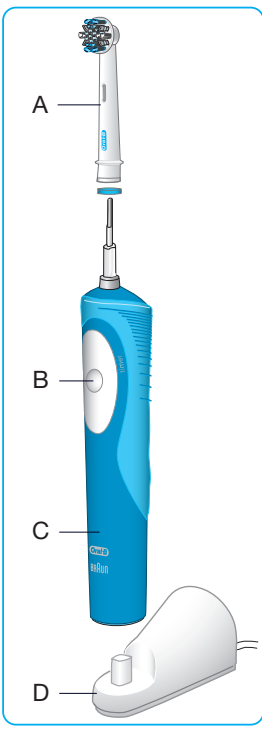IMPORTANT SAFEGUARDS

Your Oral-B Vitality toothbrush has been developed with leading dental professionals to effectively remove plaque from all areas of your teeth.

Read all instructions before using.

DANGER

In order to reduce the risk of electrocution:

1. Do not place or store the charger where it can fall or be pulled into a tub or sink.
2. Do not place the charger in water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.
4. Do not use while bathing.

WARNING

To reduce the risk of burns, electrocution, fire or injury:

1. Close supervision is necessary when this product is used by, on, or near children or invalids.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Keep cord away from heated surfaces.
4. Never drop the appliance or insert any object into any opening of the appliance.
5. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged or if it has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.

IMPORTANT

1. Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Center. A damaged or non-functioning unit should no longer be used.
2. If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible
3. Do not touch the power plug with wet hands. This can cause an electric shock.
4. When unplugging, always hold the power plug instead of the cord.
5. Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Service center. Do not disassemble the product

except when disposing of the battery.

6. When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
7. This product is not intended for use by children under age three.

CAUTION

If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

SAVE THESE INSTRUCTIONS

Description

- A Brush head (depending on model) – interchangeable
- B On/off switch
- C Handle
- D Charging unit

Specifications

For voltage specifications please refer to the bottom of the charging unit.

Connecting and charging

Your toothbrush has a water-proof handle, is electrically safe and designed for use in the bathroom.

- Plug the charging unit (D) into an electrical outlet. Place the handle (C) on the charging unit.
- A full charge takes 16 hours and enables approximately one week of regular brushing (twice a day, 2 minutes).
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible. However for environmental reasons, Oral-B recommends to unplug the charging unit until the next recharge is required.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

Using the toothbrush

- Wet brush head and apply toothpaste. You can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on (B).
- Slowly guide the brush head from tooth to tooth, spending a few seconds on each tooth surface.
- Brush the gums as well as the teeth, first the outsides, then the insides and finally the chewing surfaces (1–5). Brush at least 2 minutes to thoroughly remove plaque. Do not press too hard or scrub.
- Turn off the handle by pressing the on/off switch.

Timer (depending on model)

An incorporated timer memorises the elapsed brushing time. A short stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

Individual brush heads and accessories

Oral-B offers you a variety of brush heads which are available in stores that sell Oral-B products. Following brush heads are interchangeable and can be used on your Oral-B Vitality handle.

Oral-B FlossAction brush head

The Oral-B FlossAction brush head is recommended for everyday cleaning. It

features micropulse bristles that get in hard-to-reach areas

between your teeth, gently

removing hidden plaque.

Oral-B Precision Clean

brush head

The Oral-B Precision Clean brush head provides effective and gentle cleaning

action for both teeth and gums.

The brush head cups each

tooth with longer and denser

bristles and reaches deep

between teeth at the same

time.

Oral-B Sensitive brush head

The Oral-B Sensitive brush head provides thorough cleaning

for sensitive teeth

and gums. It has an extra soft

bristle texture that is gentle on

teeth and gums.

Oral-B ProWhite brush head

For occasional or everyday

polishing, please

use the Oral-B

ProWhite brush

head. It features a dentist-

inspired polishing cup to

remove surface stains from

coffee, tea and tobacco for

effective brightening. We do

not recommend using the

Oral-B ProWhite brush head

with braces.

Oral-B Dual Clean brush head

The Oral-B Dual Clean brush head

offers two types of cleaning

motion in one. It combines

the oscillating and

pulsating motion

with an additional side-to-side

motion. The oscillation motion

stimulates gums while remov-

ing plaque and bacteria. Multi-

angled bristles clean between

teeth and along the gumline.

You can use this brush head

either by moving from tooth

to tooth or apply the same

brushing technique as you do

with any manual toothbrush.

All brush heads feature light

blue INDICATOR™ bristles to

help you monitor brush head

replacement need. With

thorough brushing, twice a day

for two minutes, the blue color

will fade halfway down

approximately in a 3-month

period (l.b). If the bristles splay

before the colour recedes (l.c),

you are exerting too much

pressure on teeth and gums.

In the first few days of using

any electric toothbrush, your

gums may bleed slightly. In

general bleeding should stop

after a few days. Should it

persist after 2 weeks, please

consult your dentist or dental

hygienist.

Cleaning


After use, rinse the brush head for several seconds under running water with the handle switched on (II). Then switch off the handle and take off the brush head. Rinse both parts separately under running water (III) and wipe them dry. From time to time, also clean the charging unit with a damp cloth (IV).

Environmental notice

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Service Center or at appropriate collection points provided in your country.

Battery removal at the end of the product's useful life

Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

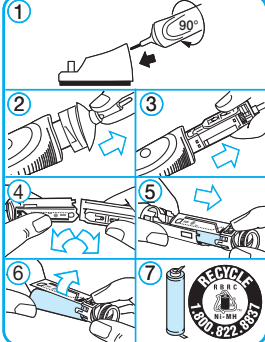Caution! Opening the handle will destroy the appliance and invalidate the warranty.

Subject to change without notice.

Satisfaction Guaranteed!

Take the challenge! Experience the Oral-B Difference! We guarantee that you will be completely satisfied.

Try any Oral-B Vitality rechargeable toothbrush for 30 days, starting from the day of purchase. If you are not 100% satisfied, receive a full refund. Guaranteed. For a full refund, please return your handle, refills and charging unit with the original cash register receipt to the address below postmarked within 30 days of purchase.

U.S. only:

Attention: 30-Day Challenge
Oral-B Braun Division
The Procter & Gamble
Company

Dept. R

PO Box 900022

El Paso, TX 88590-0022

Canada only:

Attention: 30-Day Challenge
Braun Canada a division of
Procter & Gamble Inc.

P.O. Box #64

Brampton, ON L6V 2K7

The shipment cost to return the product is the responsibility of the consumer.

Terms and conditions apply, for further details, please contact our customer service center.

FOR USA ONLY

Limited 2-Year Warranty

Oral-B warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Oral-B will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Oral-B Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Oral-B. This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL ORAL-B BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

How to obtain service in the USA

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Service Center of your choice or,
- B. Ship the product to the authorized Oral-B

Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Oral-B replacement parts and accessories:

Contact an authorized Oral-B Service Center. Please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To order original Oral-B replacement parts and accessories (including replacement brush heads) by phone:

Please call 1-800-566-7252.

To speak to an Oral-B Consumer Service Representative:

Please call 1-800-566-7252.

FOR CANADA ONLY

Limited 2-Year Warranty

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service center, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages. The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

How to obtain service in Canada

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Oral-B replacement parts & accessories:

Contact an authorized Oral-B Service Center. Please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To speak to an Oral-B Service Representative:

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Braun GmbH